L A U N C H

FGM

presents:

Elóra

A CONTEMPORARY COMPOSITE

ELÓRA APS TECHNICAL GUIDE

Our goal is to ensure that you achieve the best results when using **Elóra APS** in your restorative procedures and make the most of all the benefits that the composite offers.

To this end, we have created this technical guide, in which you will find a list of clinical indications, its shade system, and a step-by-step guide with restorative suggestions for the most common procedures, as well as clinical cases and product usage instructions.

Discover the contemporary nature of **Elóra APS** and its various possibilities.


FGM
Elóra
APS
Resina composta fotopolimerizável
Light curing composite resin
Resina compuesta fotopolimerizable
FGM
Elóra
APS



Photo: Dr. Paulo Gabriel Warmeling



CONTENT

Introduction 06

Shade System 08

Indications 10

Restorative Possibilities 11

Technology for High-Esthetic Results 12

Excellent Handling 13

APS Technology 13

Mechanical Properties 14

CLINICAL CASES 18

Esthetic contouring in a diastema case associated with enamel defect 18

Diastema closure using a simplified shade system 24

Class II Restoration with Elóra APS composite 32

Composite veneers for bleached teeth with Elóra APS 36

Restoring dental esthetics with diastema closure 42

Commercial Presentation 51



Photo: Dr. Ariel Braga

A contemporary composite.
One concept, many possibilities.

Prof. Dr. Maria Carolina Erhardt
Associate Professor of operative dentistry at the Faculty of Dentistry, Federal University of Rio Grande do Sul

Prof. Dr. Thais Tomé Feldens
Associate Professor of operative dentistry at the Faculty of Dentistry, Federal University of Rio Grande do Sul

The new **Elóra APS** composite brings an innovative proposal with its simplified shade system and advanced composition. This nanohybrid composite, with a differentiated blend of fillers and APS curing technology, combines high resistance, excellent polishing, and durability, precisely meeting the most prevalent clinical demands.

The major advantage of this simplified composite system is the ease of shade selection. With only 12 shades, covering the 19 tones of the VITA shade guide, the choice becomes more practical and intuitive. The reduced palette, consisting of shades with intermediate translucency, allows predictable and natural aesthetic results, avoiding the need for multiple layers. In addition to reducing clinical time, there is an increase in the accuracy of the final restoration shade.

The versatility of the **Elóra APS** system stands out by allowing a single material to be applied in different types of restorations, simplifying the workflow and reducing the need for large material stocks. This brings significant benefits both for office management and cost optimization, directly impacting clinical planning and treatment efficiency.

Another aspect worth highlighting is the APS curing technology, which extends working time and provides optimized shade stability over time. This is particularly advantageous in situations that require precise handling and long-lasting aesthetic finishing. Furthermore, the dentin shades of the **Elóra APS** system prove effective in masking darkened substrates, offering excellent aesthetic results even in challenging conditions.

The nanohybrid composition of **Elóra APS** also provides high wear resistance, an essential characteristic for composites that aim to last and maintain their appearance during prolonged clinical function. At the same time, the material is easy to handle, making the application and finishing process more comfortable for both the clinician and the patient.

Elóra APS represents an advancement in Restorative Dentistry, combining practicality, aesthetic quality, and functionality in a system that simplifies and optimizes the clinical routine.



SHADE SYSTEM

With **Elóra APS**, it is possible to achieve high-level aesthetics in a simple and effective way. Its specially developed composition features a modern and optimized shade system that delivers impressive results, with translucency and optical properties balanced with natural teeth.

In either the layering or the simplified technique, **Elóra APS** meets a wide range of clinical applications.

SYSTEM WITH DIFFERENT LEVELS OF TRANSLUCENCY

BODY

XBL BL B1 A1 A2 A3 A3.5

DENTIN

DB1 DA1 DA3

EFFECT ENAMEL

TRANS WE









CATEGORY	DENTIN	BODY	EFFECT ENAMEL	
SHADES	DB1, DA1 E DA3	A1, A2, A3, A3.5, B1, BL, XBL	TRANS	WE
DESCRIPTION	Composites with lower translucency, indicated for use as the first layer in dentin reproduction.	Composites with medium translucency that can be used as an intermediate layer or as a single layer, reproducing both enamel and dentin.	The most translucent effect composite of the **Elóra APS** system, used to reproduce effects in specific areas such as incisal and proximal regions.	The achromatic enamel composite can be used, depending on the type of restoration, as the palatal and buccal enamel, or only as the buccal enamel (for example, in veneer cases).
CURING TIME	30 Seconds	20 Seconds	20 Seconds	



Photo: Dr. Paulo Gabriel Warmling


INDICATIONS

Wide range of clinical indications for anterior and posterior teeth

- ✓ Direct restorations of all classes;
- ✓ Esthetic/diastemata recontouring;
- ✓ Core build-ups;
- ✓ Indirect and semi-direct restorations;
- ✓ Porcelain/composite repairs.

RESTORATIVE POSSIBILITIES

Its presentation is suitable for the main cases encountered in daily clinical practice, offering high mechanical strength for posterior teeth and meeting the highest aesthetic demands of anterior teeth.

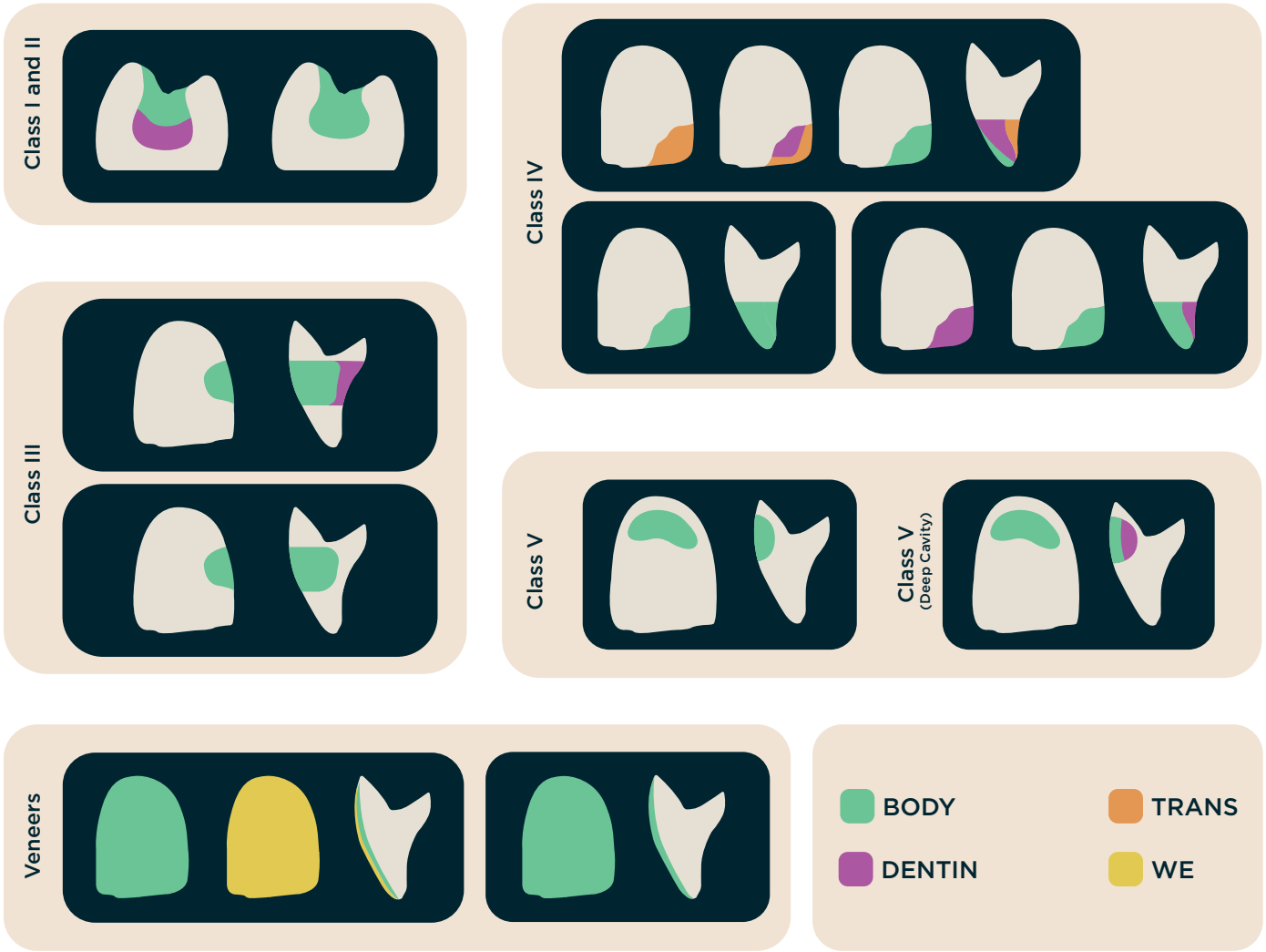

TECHNOLOGY FOR HIGH-AESTHETIC RESULTS


Elóra APS is a radiopaque composite, featuring a high-performance monomer composition and a new combination of fillers, with a wide distribution of particle sizes for greater packing density, forming the **Optimized Filler Distribution technology.**

4 TYPES OF **FILLERS**

- Barium Glass;
- Spherical Zirconium Silicate;
- Spherical Synthetic Silica;
- Amorphous Synthetic Silica.

DISTRIBUTED IN **6** DIFFERENT **PARTICLE SIZES**

Promoting optimized filler distribution and packing from 0.02 to 0.7 microns within the polymer.

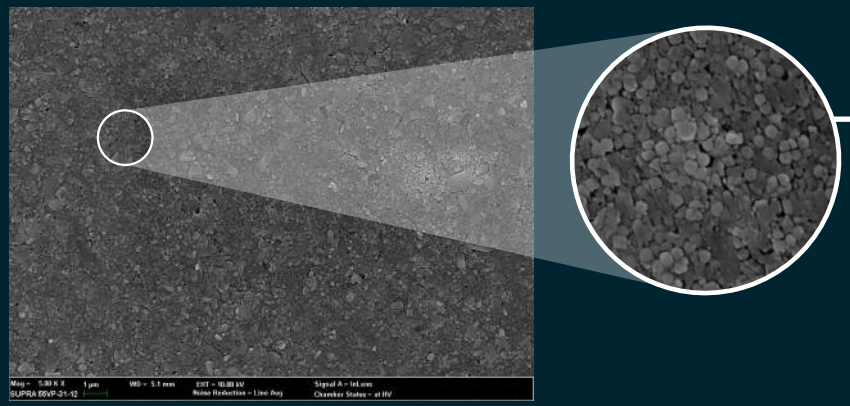

Source: Internal data FGM Dental Group.

The combination of

POLYGONAL FILLERS — SPHERICAL FILLERS

30% **70%**


Result in better packing and condensation, providing the product with:

OPTIMIZED FILLER DISTRIBUTION

SPHEROIDAL ZIRCONIA SILICATE

FILLER CONTENT **+76%**

THE BEST MECHANICAL PROPERTIES AND HIGH GLOSS.

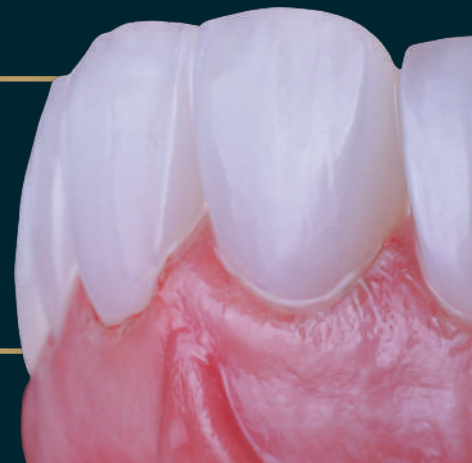

EXCELLENT HANDLING

The composition balance of Elóra APS composite optimizes the product's viscosity while sculpting, resulting in an easy and pleasant handling consistency that does not stick to the spatula and adapts easily to the surface.

These characteristics make polishing the composite simple, achieving excellent surface gloss.

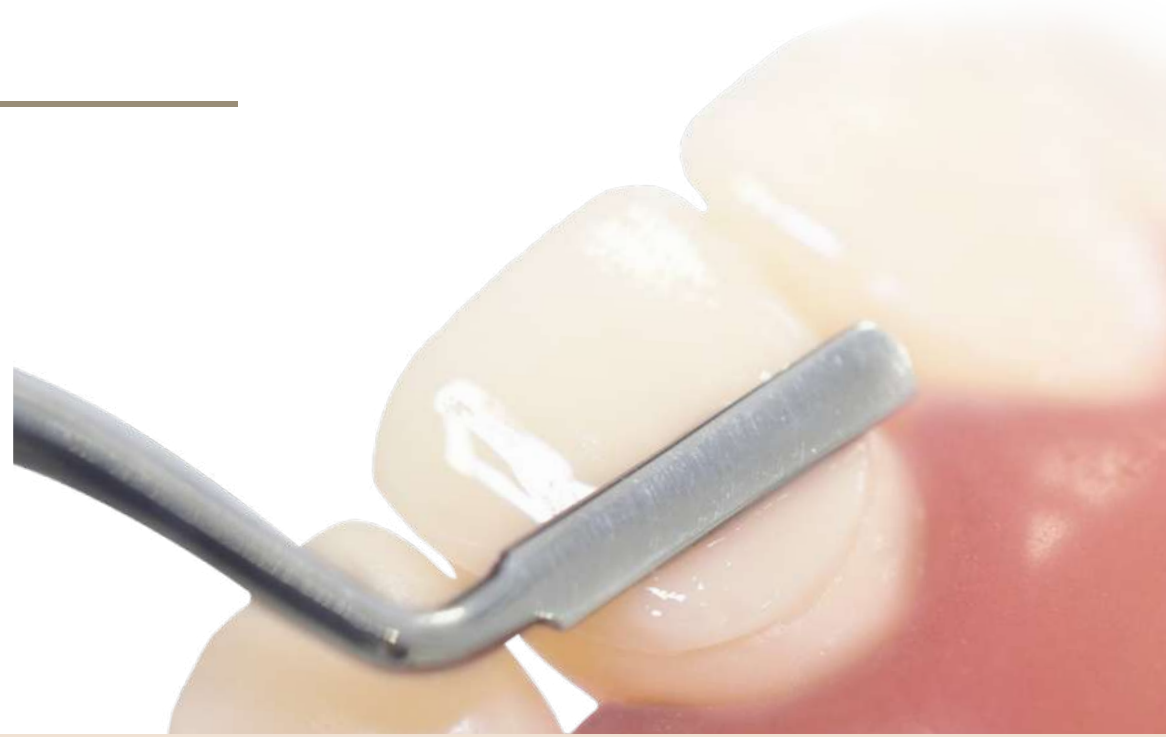


Photo: Dr. Orlando Reginatto



APS TECHNOLOGY


FGM Exclusive

- Provides more time for modeling and sculpting the composite;
- Shade predictability, with no significant variation after curing;
- Excellent mechanical properties.

Technology APS ADVANCED POLYMERIZATION SYSTEM

MECHANICAL PROPERTIES

Roughness (% Ra)

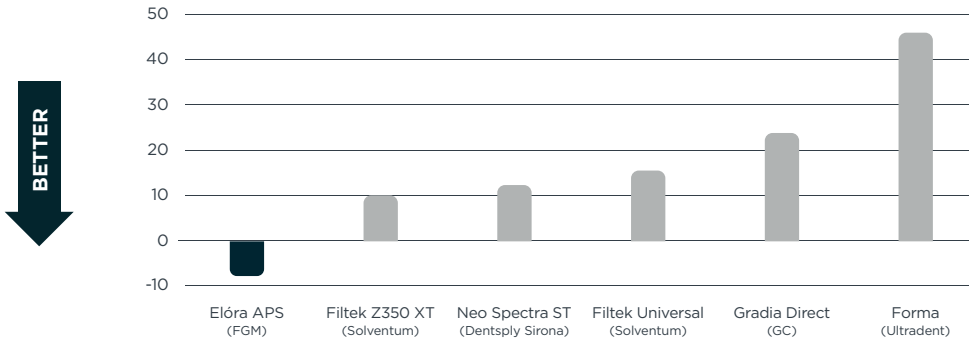


Graph 1. Percentage of roughness when comparing values before and after simulated brushing (n=10 per experimental condition).

Source: Moreira PHA, Cochinski GD, Sovinski G, Reis A, Bauer J, Loguercio A. State University of Ponta Grossa (UEPG) and Federal University of Maranhão (UFMA), 2024.


Elóra APS was the **ONLY** to not show an increase in surface roughness after simulated brushing, when compared to competing composites.

Knoop Microhardness (KHN)

Knoop Microhardness (Surface) - 2A

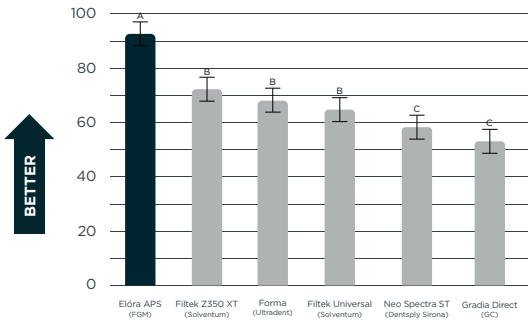

Knoop Microhardness (Bottom) - 2B

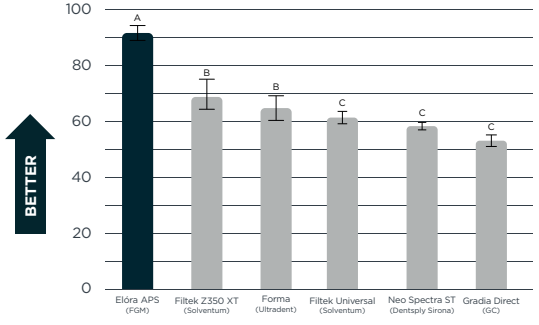

Graphs 2A and 2B. Mean and standard deviation of Knoop microhardness (KHN) of the composites evaluated (n=10 per experimental condition).

(*) Identical letters indicate statistical similarity (2-way ANOVA and Tukey's test; $p<0.05$)
Source: Moreira PHA, Cochinski GD, Sovinski G, Reis A, Bauer J, Loguercio A. State University of Ponta Grossa (UEPG) and Federal University of Maranhão (UFMA), 2024.


Elóra APS showed the **HIGHEST** Knoop microhardness values (surface and bottom), with statistically **SUPERIOR** results when compared to the evaluated competing composites.



Benefits in clinical practice

In practice, these characteristics result in greater resistance of the restorative material to wear and abrasion, contributing to the longevity of the restoration. Furthermore, composites with these properties, such as Elóra APS, tend to better preserve their gloss and smoothness over time, ensuring more durable aesthetics.

Flexural Strength (MPa)

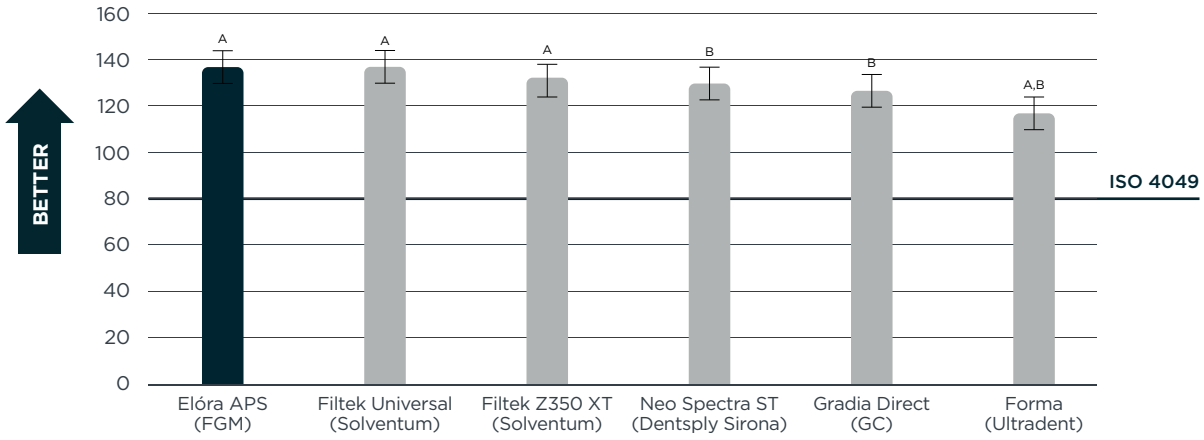

Graph 3. Mean and standard deviation of flexural strength (MPa) of the composites evaluated (n=10 per experimental condition).

(*) Identical letters indicate statistical similarity (2-way ANOVA and Tukey's test; $p<0.05$)
Source: Moreira PHA, Cochinski GD, Sovinski G, Reis A, Bauer J, Loguercio A. State University of Ponta Grossa (UEPG) and Federal University of Maranhão (UFMA), 2024.

Flexural Strength is a mechanical property that allows the evaluation of the composites' resistance to fracture. **Elóra APS** demonstrated a high flexural strength value, well above the established standard for composites, positioning it among the best composites on the market.




Benefits in clinical practice

In practice, it is essential to select materials with adequate mechanical strength, such as Elóra APS, especially for teeth subjected to high masticatory loads, such as occlusal surfaces of posterior teeth and incisal edges of anterior teeth. In this regard, flexural strength reflects the quality of the composite formed after curing, resulting from the interaction between the polymer matrix, the types of fillers used, and the efficiency of photoinitiators in converting monomers into polymers.

Dentin Bond Strength (MPa):

Association of Elóra APS composite with different adhesives

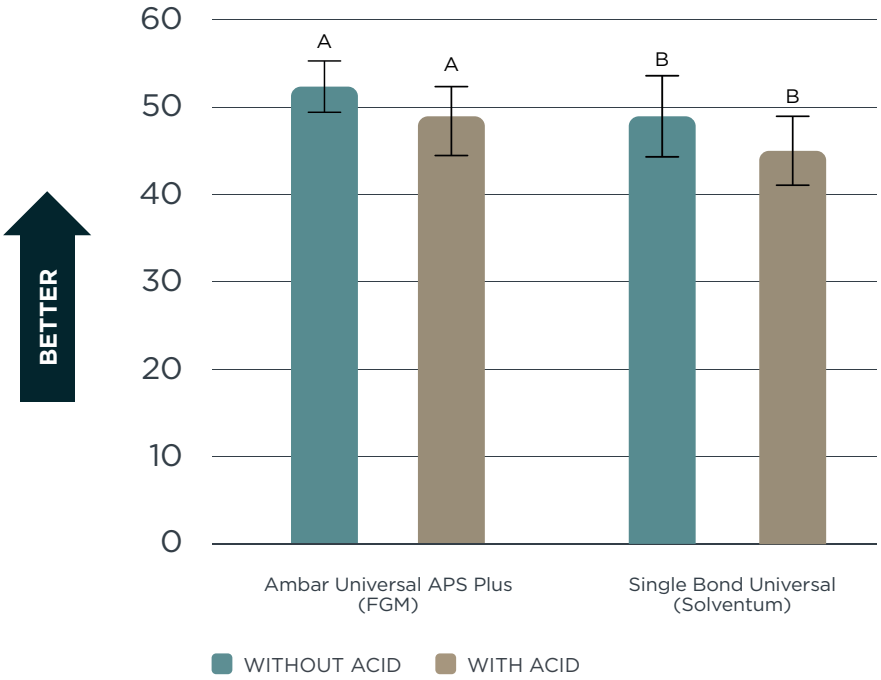

Graph 4. Mean and standard deviation of dentin bond strength (MPa) of Elóra APS composite to dentin mediated by different adhesive systems (n = 5 per experimental condition) (*).

(*) Identical letters indicate statistical similarity (1-way ANOVA and Tukey's test; $p < 0.05$).
Source: Sovinski G, Moreira PHA, Barbosa LMB, Cochinski GD, Reis A, Bauer J, Loguercio A. State University of Ponta Grossa (UEPG) and Federal University of Maranhão (UFMA), 2024.

The association between **Elóra APS** composite and the Ambar Universal APS Plus adhesive system (FGM) showed superior adhesion to dentin (with and without phosphoric acid) when compared to **Elóra APS** composite associated with the Single Bond Universal adhesive system (Solventum).




Benefits in Clinical Practice

In practice, high dentin bond strength values indicate the ability of the adhesive material (such as an adhesive system and composite) to form a lasting and stable bond with the dentin substrate over time, resulting in reduced marginal infiltration.

Dentin Bond Strength (MPa):

Ambar Universal APS Plus associated with different composites

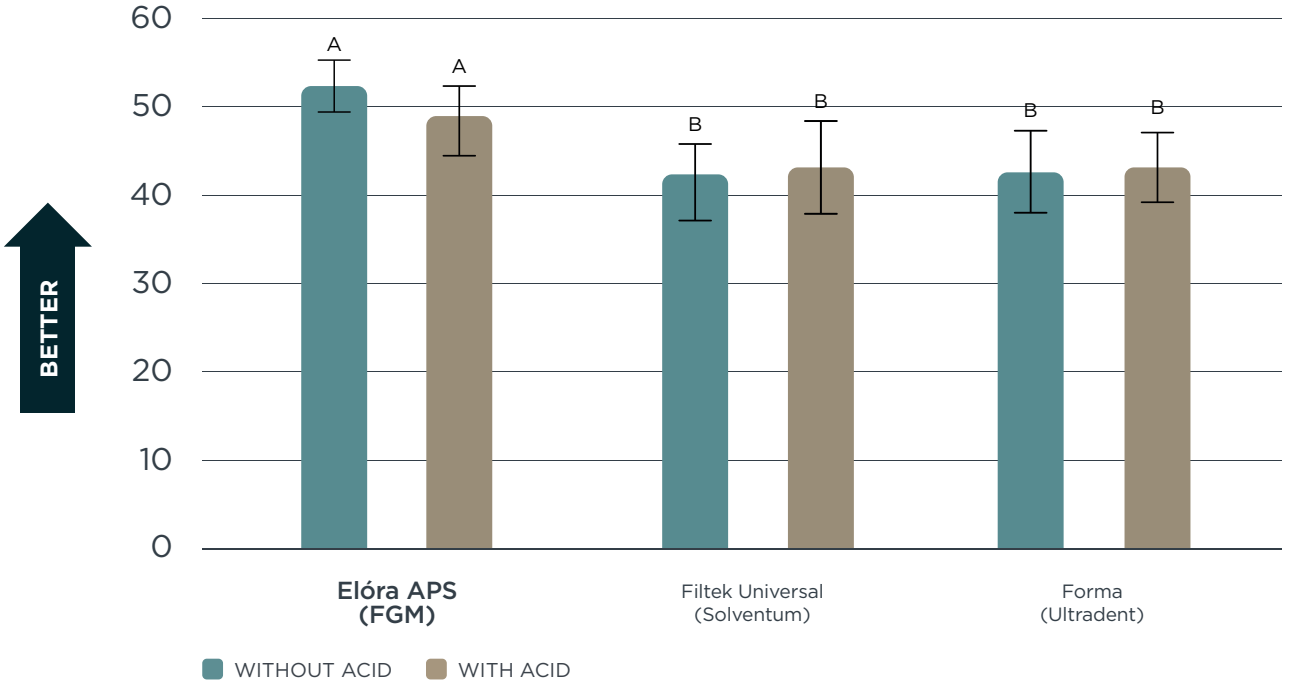

Graph 5. Mean and standard deviation of dentin bond strength (MPa) of the Ambar Universal APS Plus adhesive associated with different composites, with and without phosphoric acid (n = 5 per experimental condition) (*).

(*) Identical letters indicate statistical similarity (1-way ANOVA and Tukey's test; $p < 0.05$).
Source: Sovinski G, Moreira PHA, Barbosa LMB, Cochinski GD, Reis A, Bauer J, Loguercio A. State University of Ponta Grossa (UEPG) and Federal University of Maranhão (UFMA), 2024.

The association between **Elóra APS** composite and the Ambar Universal APS Plus adhesive system showed statistically **SUPERIOR** dentin bond strength results when compared to the evaluated competing composites.

For optimal bond strength results and restoration longevity, it is recommended to use products from the same FGM line.

CLINICAL CASE 01

Author: Prof. Paulo Gabriel Warmling | Female Patient, 25 years old.



01. Initial appearence of the smile, showing the presence of diastemas and enamel defect on tooth 11.

Aesthetic contouring in a diastema case associated with enamel defect.

Main complaint: Aesthetic dissatisfaction due to stain on anterior tooth and tooth shape.

Case Report: The patient sought treatment due to dissatisfaction with the spaces between the anterior teeth and a small enamel defect on the central incisor, which compromised the aesthetics of the smile.

Initially, the operative field was isolated with a rubber dam to ensure a dry and controlled environment, acid etching was then performed on the enamel using **Condac 37**, followed by active application of the **Ambar Universal APS Plus** adhesive, which was light-cured for 10 seconds with the **Quazar** device.


After the adhesive protocol, the stratification stage began. For this case, only two shades of **Elóra APS** composite were used, which allowed the masking and excellent coverage of the stain and enamel defect on tooth 11 in a conservative manner. Additionally, aesthetic contouring was performed to close the diastemas on teeth 12, 21, and 22 using the same shades. After stratification, the finishing protocol was carried out using the **Diamond Pro** abrasive discs, followed by final polishing with the **Diamond Flex** felt disk and **Diamond Excel Ultrafine** polishing paste.



02




03

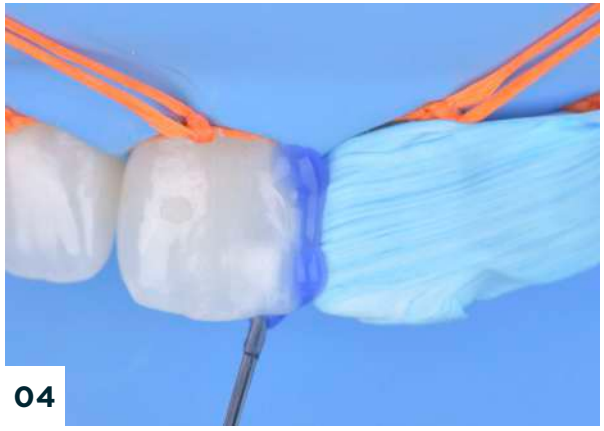


04



05A

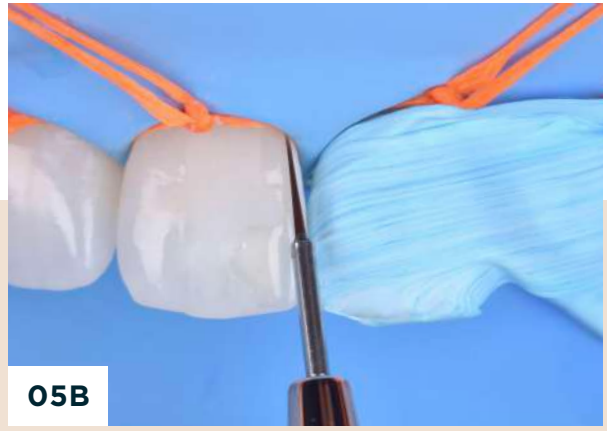

05B

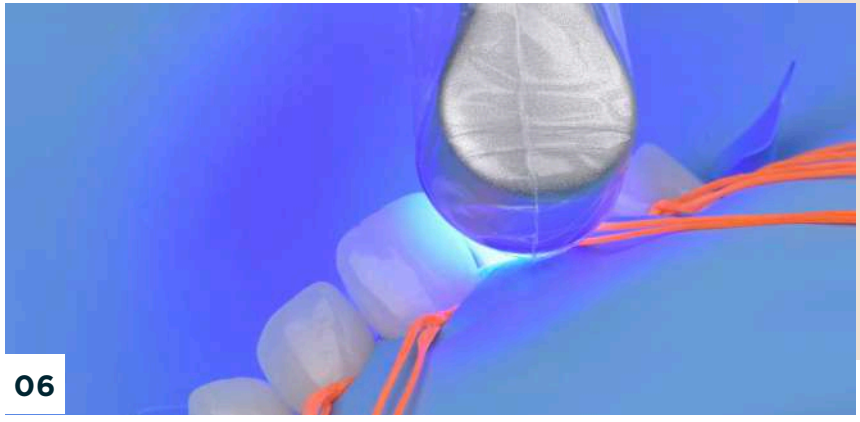


06

02. Intraoral appearence.
03. Absolute isolation.
04. Enamel etching on the mesial of tooth 11 with **Condac 37.**
05A and 05B. Application of **Ambar Universal APS Plus** adhesive on enamel.
06. Light-curing of the adhesive with the **Quazar** device.



07

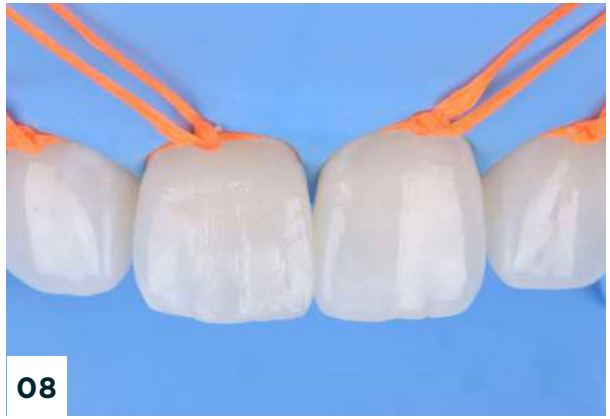

08



09



10A



10B

07. Interproximal adaptation of the polyester strip and composite increment at the cervical of tooth 11 with **Elóra APS**, in shade B1.
08. Appearence after composite application and light-curing.
09. Application of mesio-vestibular increment on tooth 21 with **Elóra APS** composite, in shade B1.
10A and 10B. Application of mesial increment on tooth 22 with **Elóra APS** composite, in shade BL, and appearence after contouring.




11A



11B



12

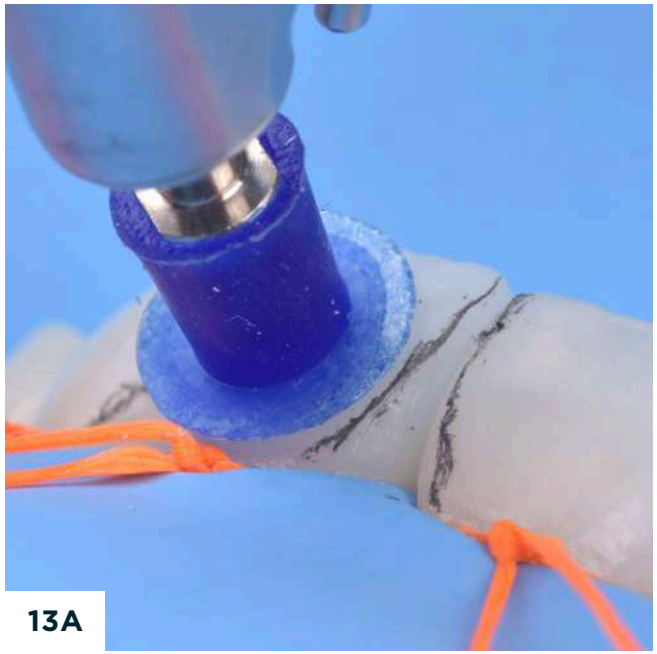

13A



13B

11A and 11B. Application of **Elóra APS** composite, in shade BL on the distal of tooth 11 and appearence after aesthetic contouring.
12. Application of **Elóra APS** composite, in shade BL on the mesio-vestibular surface of tooth 12.
13A and 13B. Finishing with **Diamond Pro** sanding discs, final polishing with **Diamond Flex** felt disc, and **Diamond Excel Ultrafine** polishing paste.



14A



14B



14C

14A, 14B and 14C. Final appearance after aesthetic contouring – front view and lateral views.



15

15. Final appearance.



Didactic diagram of stratification.

- Elóra APS, in shade B1.
- Elóra APS, in shade BL.

CLINICAL CASE 02

Author: Dr. Ariel Braga | Female Patient, 26 years old.



01

01. Initial appearance of the smile, frontal view showing diastemas between the anterior teeth.

Diastema closure using a simplified shade system

Main complaint: Aesthetic dissatisfaction with spaces between the anterior teeth.

Case report: The patient reported aesthetic dissatisfaction regarding the spaces between the anterior teeth. After the anamnesis and clinical examination, it was found that the gum health was good, and there were no pathologies, making aesthetic restorations possible according to the patient's wishes.

Initially, a prophylaxis was performed, followed by relative isolation and the insertion of retraction cord for moisture control and gingival retraction. After testing for shade and translucency, it was decided to use the **Elóra APS** composite, in shade A2. Its balanced translucency allowed the diastema closures to be effectively and naturally completed, without any grayish appearance.

After the diastema closures on teeth 11, 12, 21, and 22 were completed, the finishing and polishing protocol was carried out using the **Diamond Pro** series of polishing discs and felt discs, along with **Diamond Flex** polishing paste and **Diamond Excel Ultrafine** polishing paste. The result observed was of extreme naturalness, with excellent gloss and surface smoothness, meeting the aesthetic requirements and enabling a simplified stratification protocol.



02



03




04



05

02. Initial appearance of the smile, lateral view.
03. Intraoral appearance.
04. Relative isolation, and shade and opacity test with **Elóra APS** composite, body A1 and A2.
05. Insertion of retraction cord for gingival retraction.

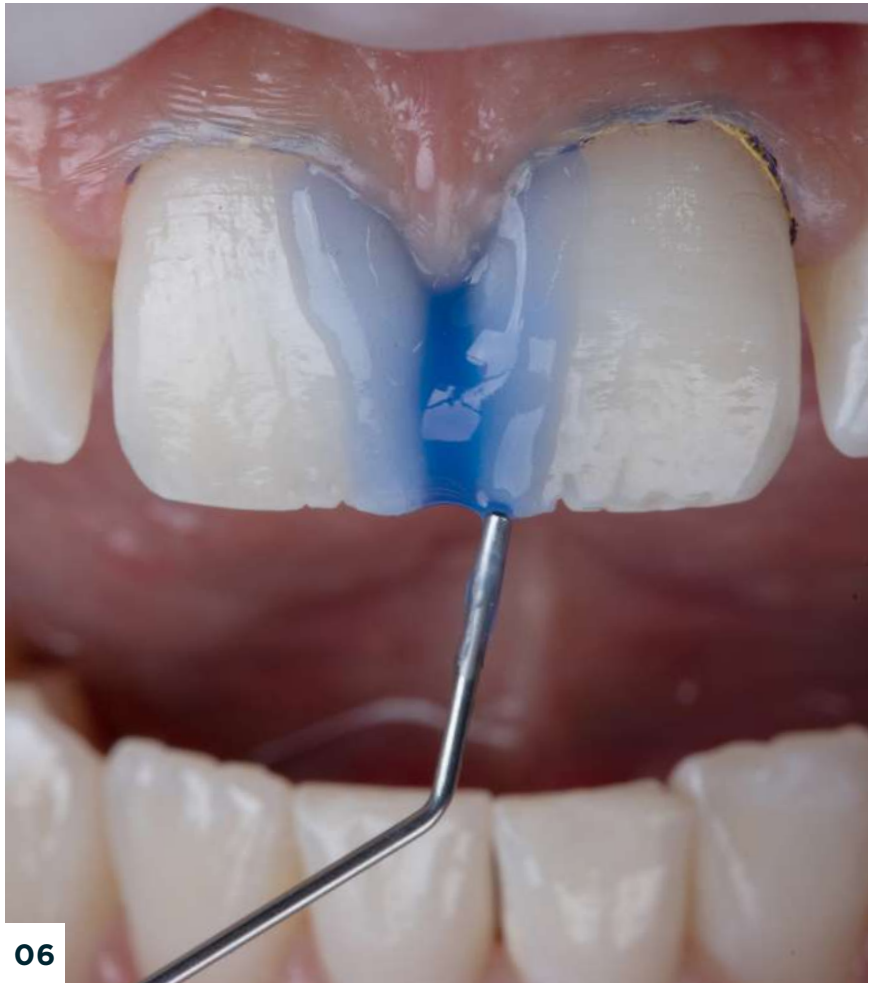

06




07



08


06. Acid etching on enamel with **Condac 37.**
07. Appearance after enamel etching.
08. Application of **Ambar Universal APS** adhesive.

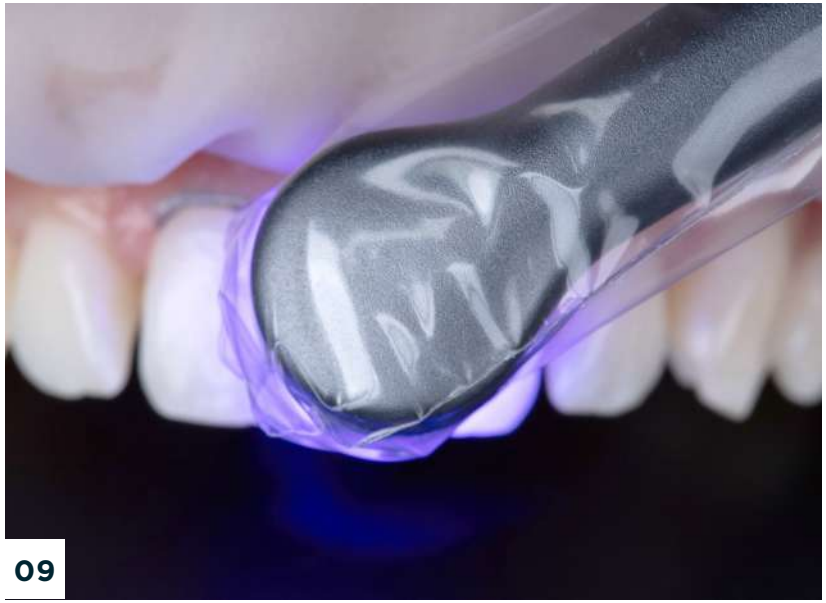

09



10

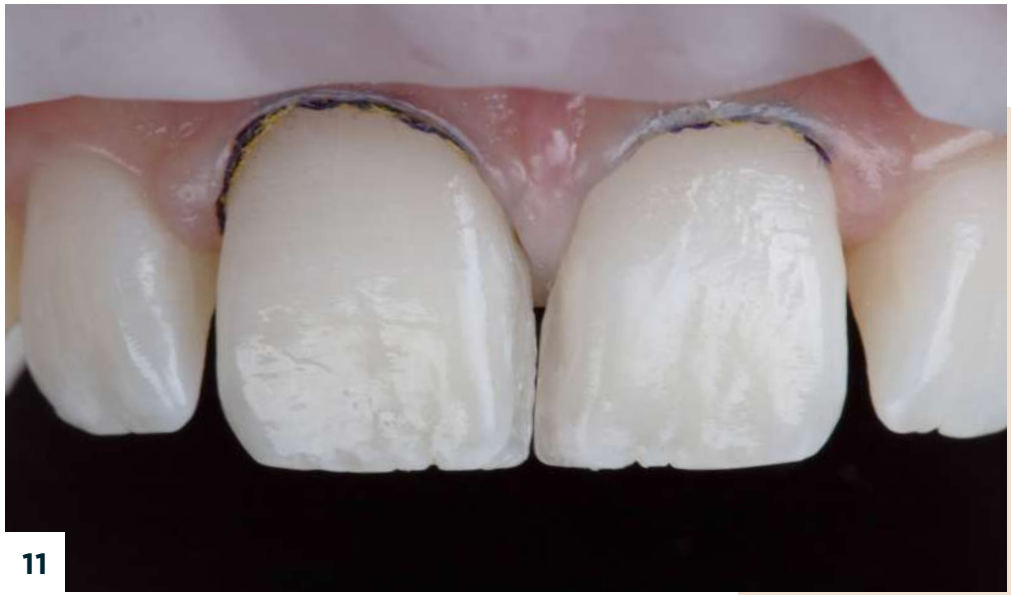

11

09. Light-curing of the adhesive with **Quazar** for 10 seconds.
10. Construction of the palatal wall to close the diastema on tooth 11 with **Elóra APS** composite, in shade A2.
11. Construction of the palatal wall to close the diastema on tooth 21 with **Elóra APS** composite, in shade A2.



12



13



14

12. Application of the vestibular layer with **Elóra APS** composite, in shade A2.
13. Appearance after light-curing.
14. Finishing with the **Diamond Pro** polishing disc sequence.



15



16



17

15. Final polishing with **Diamond Flex** felt discs and **Diamond Excel Ultrafine** polishing paste.
16. Construction of the palatal layer for diastema closure on tooth 12 with **Elóra APS** composite, in shade A2.
17. Application of the vestibular layer on tooth 12 with **Elóra APS** composite, in shade A2.



18

18. Intraoral appearance.
19. Final smile appearance, lateral view.



19



20

20. Final appearance.



Didactic diagram of stratification.

Elóra APS, in shade A2.

CLINICAL CASE 03

Author: Prof. Dr. Renato Voss Rosa | Male Patient, 37 years old.



01

01. Initial appearance of the previous restoration.

Class II Restoration with Elóra APS composite.

Main complaint: Posterior tooth with unsatisfactory restoration.

Case report: The patient presented to the clinic with a poorly adapted occlusoproximal restoration in posterior tooth 16, with marginal infiltration, which had been in place for approximately 20 years.

Initially, the previous restoration was removed, followed by absolute isolation to control the operative field and ensure proper gingival retraction. Next, selective enamel etching was performed using phosphoric acid **Condac 37**, and the matrix system was positioned, ensuring proper marginal adaptation and intimate contact between the tooth and matrix.

For the restorative procedure, the material of choice was the **Elóra APS** composite, which has ideal characteristics to address clinical cases like this in a simple manner with aesthetic efficiency.

For the adhesive protocol, the **Ambar Universal APS** system was chosen, which was light-cured with the **Quazar** device. Stratification began with the construction of the marginal ridge using **Elóra APS**, in shade A2, transforming the cavity into a Class I and ensuring a good contact point. A thin layer of **Opus Bulk Fill Flow APS** composite was applied to the pulpal wall, followed by a dentin layer with **Elóra APS,** in shade DA3, and the occlusal sculpting was finished with **Elóra APS,** in shade A1. Finally, finishing and polishing were performed.



02



03



04



05

02. Appearance of the cavity after removal of the restoration and edge finishing.
03. Installation of the absolute isolation.
04. Etching with phosphoric acid **Condac 37.** Above, it is possible to observe the ability of the acid not to flow and stay in the applied area, increasing the precision of the technique.
05. Installation of the chosen matrix system.


06


07


08


09


10


11


12

06 and 07. Active application of **Ambar Universal APS** adhesive.
08. Creation of the marginal ridge with **Elóra APS** composite, in shade A2.
09. Increment of **Opus Bulk Fill Flow APS** composite on the pulpal wall.
10. Increment of dentin performed with **Elóra APS** composite, in shade DA3.
11 and 12. Creation of final increment and occlusal anatomy performed with **Elóra APS** composite, in shade A1.




13

13. Final appearance of the restoration after occlusal adjustment, finishing, and polishing.

Didactic diagram of stratification.

- Elóra APS, in shade A2
- Elóra APS, in shade A1
- Elóra APS, in shade DA3

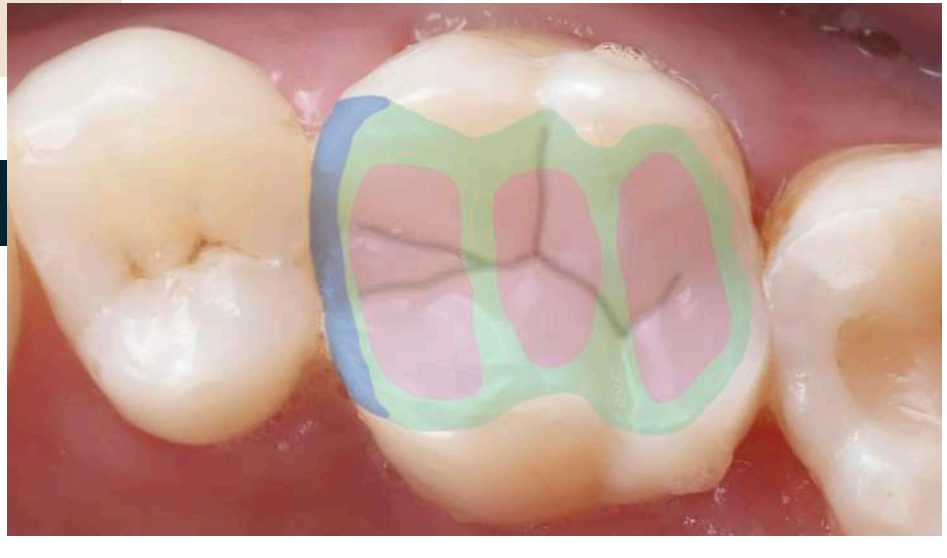

CLINICAL CASE 04

Author: Dr. Orlando Reginatto | Female Patient, 35 years old.



01

01. Initial appearance of the smile, showing disharmony in the shape of the teeth.

Composite veneers for bleached teeth with Elóra APS

Main complaint: Dissatisfaction with the harmony of the shape and shade of the anterior teeth.

Case report: The patient sought aesthetic restorative treatment because she was dissatisfied with the shape and shade of her anterior teeth. Initially, a photographic protocol was carried out for case documentation and evaluation. After a detailed planning process, it was decided in collaboration with the patient to perform composite veneers on the anterior teeth. In the initial stage, a prophylaxis was performed using pumice stone and water, followed by relative absolute isolation and shade testing of the composites to be used in all layers.

Next, gingival retraction was performed using a retractor cord, followed by a test fit of the silicone guide made with **Silic-One.** Enamel etching was performed using **Condac 37** phosphoric acid, followed by the adhesive protocol using **Ambar Universal APS Plus.**

For this case, the composite **Elóra APS** was chosen, which offers translucency and shade options ideal for the esthetics required by the clinical case. The palatal layer was made with the aid of a silicone guide using the effect resin **Elóra APS** in the shade Trans, the body layer was made with resin **Elóra APS** in the shade XBL and the final vestibular layer was used with enamel composite **Elóra APS** in the shade WE.


In the final stage, the entire finishing and polishing protocol was performed using specific drills and spirals, in addition to the polishing paste **Diamond Excel Ultrafine** and the felt disk **Diamond Flex.**


02


03


04


05

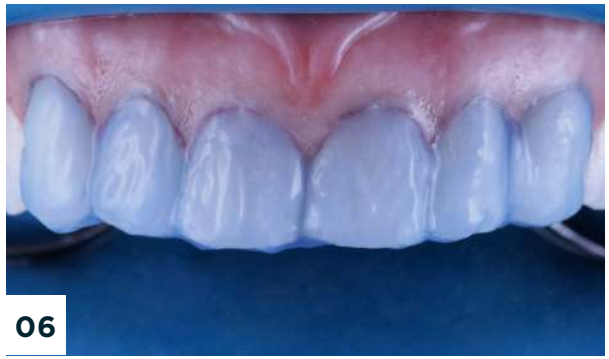
06


07


08

02. Prophylaxis with pumice stone and water.
03. After modified absolute isolation, shade testing of the **Elóra APS** composites: Trans shade for the palatal layer, shades XBL and WE for the vestibular layer.
04. Insertion of retraction cord from canine to canine.
05. Testing and marking of the palatal guide performed with addition silicone **Silic-One.**
06. Acid etching on enamel on the vestibular and palatal surfaces.
07. Active application of adhesive on enamel, with **Ambar Universal APS Plus**.
08. Light-curing with **Quazar** for 10 seconds.



09



10



11



12



13

09. Creation of the palatal layer with **Elóra APS** composite, in the Trans shade.
10. Appearance of the palatal layer after light-curing.
11 and 12. Creation of the mamelons on the incisal with **Elóra APS** composite, in the XBL shade.
13. Creation of the proximal area with the aid of polyester strip and **Elóra APS** composite, in the XBL shade.



14



15



16



17

14. Application of **Elóra APS** composite, in the Trans shade on the incisal.
15. Appearance after application of the incisal layer.
16. Application of **Elóra APS** composite, in the WE shade, on the vestibular layer.
17. Appearance and smoothness after application of the vestibular layer.









18A and 18B. Finishing protocol with diamond and multilayered drills.
19A and 19B. Polishing protocol with a sequence of rubber spirals.






20. Final polishing with **Diamond Excel Ultrafine** polishing paste and **Diamond Flex** felt disc.



21. Final appearance of the smile.

Didactic diagram of stratification.

- Elóra APS, in shade XBL
- Elóra APS, in shade TRANS
- Elóra APS, in shade WE

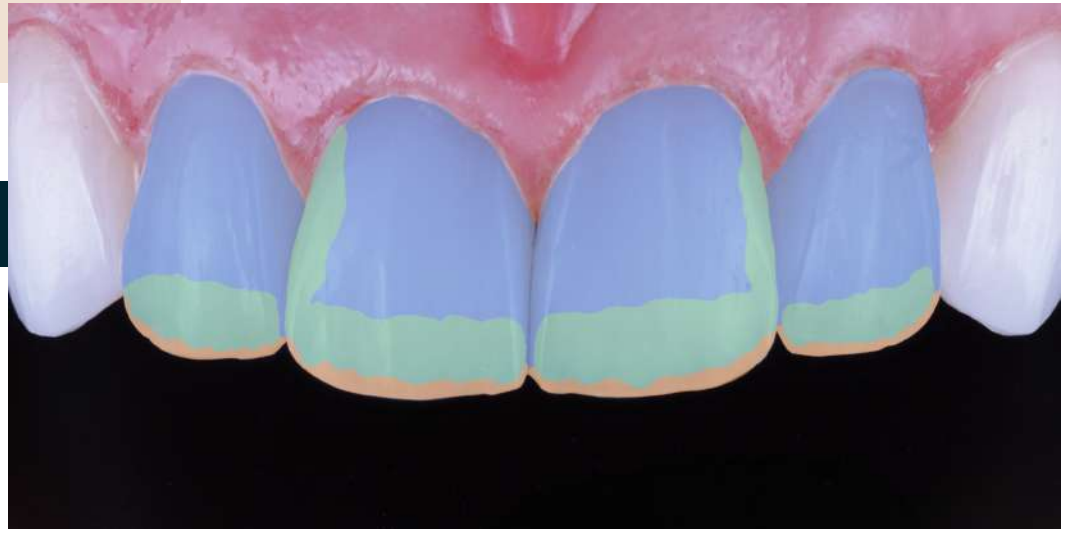

CLINICAL CASE 05

Author: Prof. Dr. Shizuma Shibata | Male Patient, 34 years old.



01

01. Initial appearance of the smile showing diastema.

Restoring dental aesthetics with diastema closure.

Main complaint: Aesthetic dissatisfaction with spaces between the teeth.

Case report: Upon clinical examination, the presence of diastemas in the upper teeth was confirmed, as reported by the patient. Initially, a shade and opacity test was performed with the composite chosen for the case, **Elóra APS,** a system that offers excellent indications for cases of small to medium diastema closure, eliminating the complexity of stratification. After prophylaxis and absolute isolation, acid etching was performed on the enamel with **Condac 37**, followed by the application of the Ambar Universal APS Plus adhesive system and photoactivation with the **Quazar** light-curing unit.


The closure of the diastemas began with tooth 21, followed by tooth 11, and subsequently teeth 13 and 23. To assist with gingival adaptation, a metallic matrix was used, along with the flowable composite **Vittra APS Unique Flow**, in combination with the composite **Elóra APS,** in shades B1 and BL. The extension of the restoration was carefully controlled, limited to the proximal portion and not exceeding the mesial crest, in order to preserve the integrity of the enamel, while respecting the natural contour of the tooth. The width of the teeth was adjusted using **Diamond Pro** sanding discs, and the restorations were completed with finishing and polishing, using rubber wheels and **Diamond Excel Ultrafine** polishing paste.




02



03

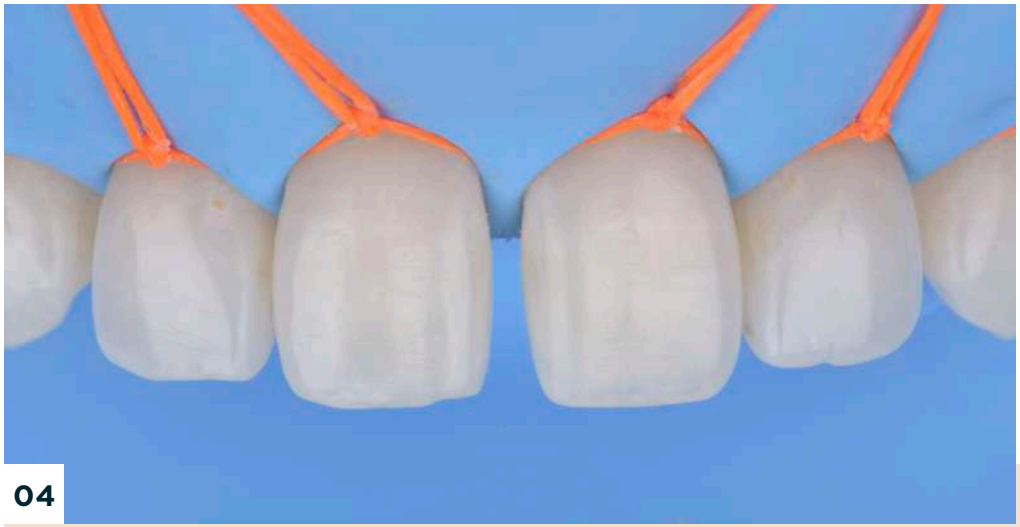


04

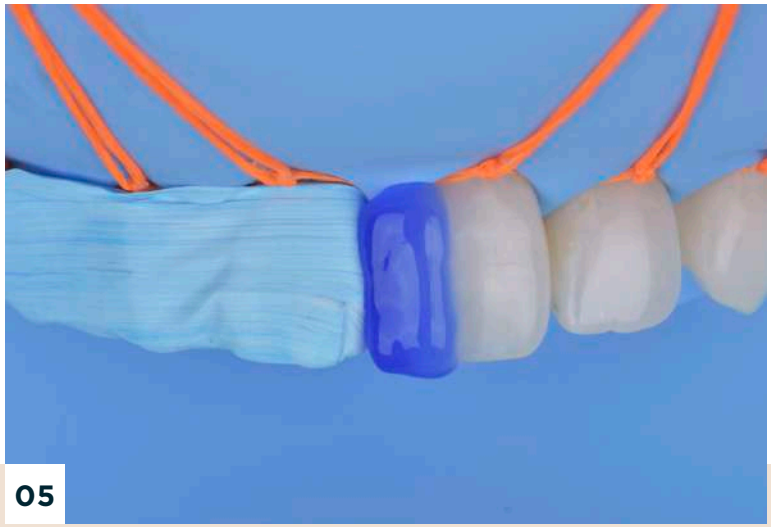


05


02. Intraoral appearance of the smile.
03. Shade and opacity test with **Elóra APS** composite, in shades B1 and BL.
04. Absolute isolation.
05. Acid etching with **Condac 37** on the mesial of tooth 21.

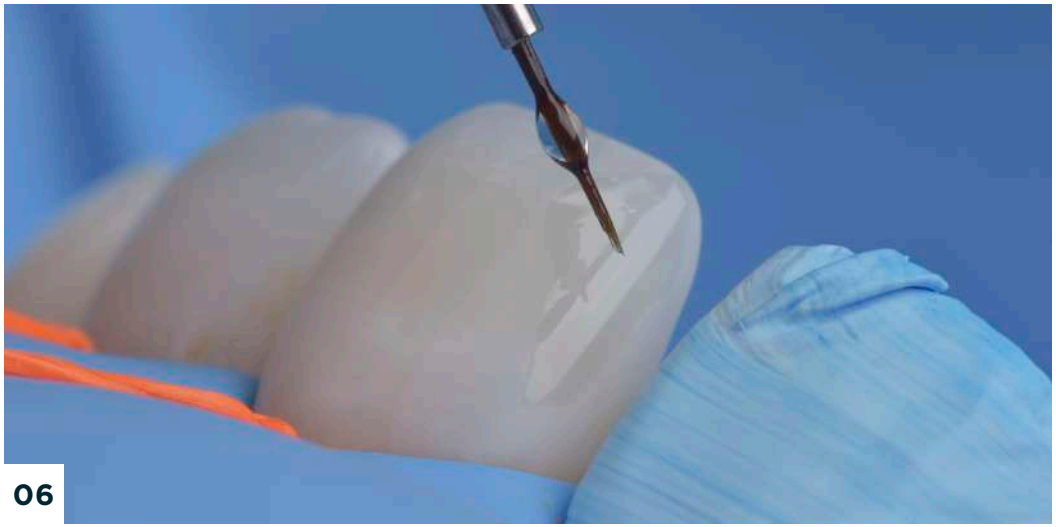


06

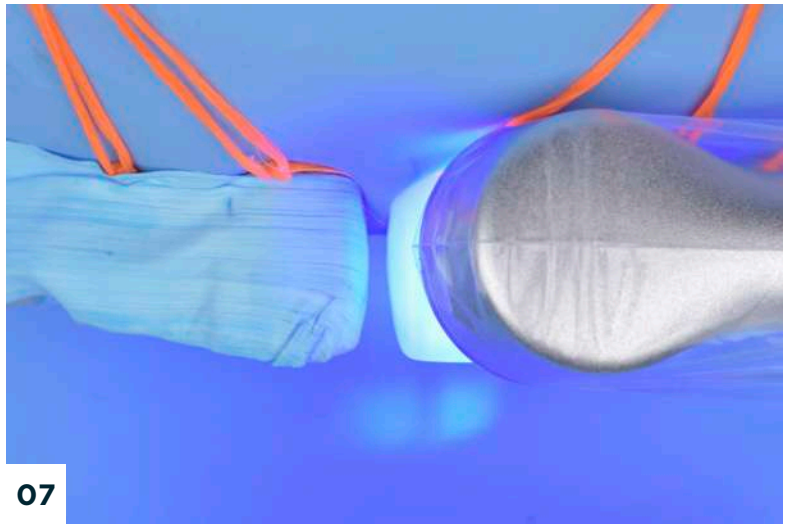


07

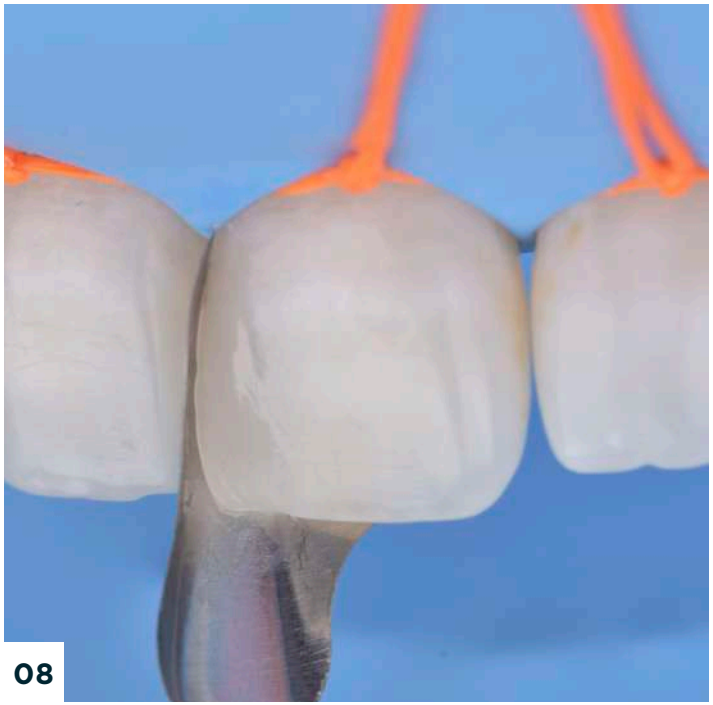

08

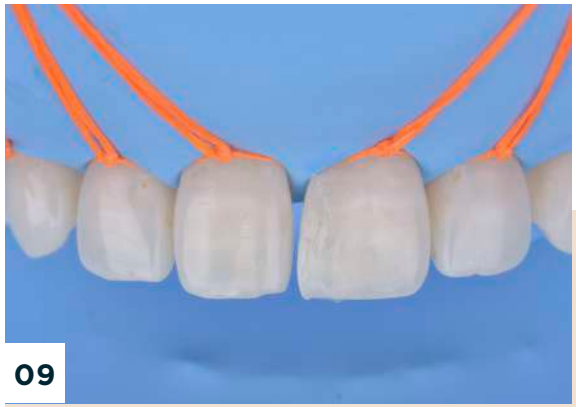


09

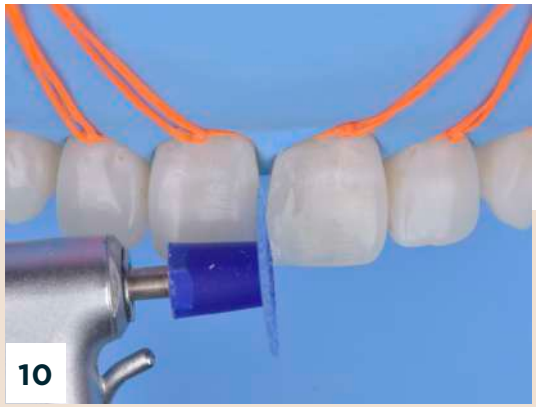


10


06. Application of **Ambar Universal APS Plus** adhesive on enamel.
07. Light-curing with **Quazar** for 10 seconds.
08. Adaptation of metallic matrix and application of **Elóra AP**S composite, in B1 shade on the cervical area and BL shade on the middle incisal third of tooth 21.
09. Appearance after matrix removal.
10. Adjustment and finishing with **Diamond Pro** sanding disc.

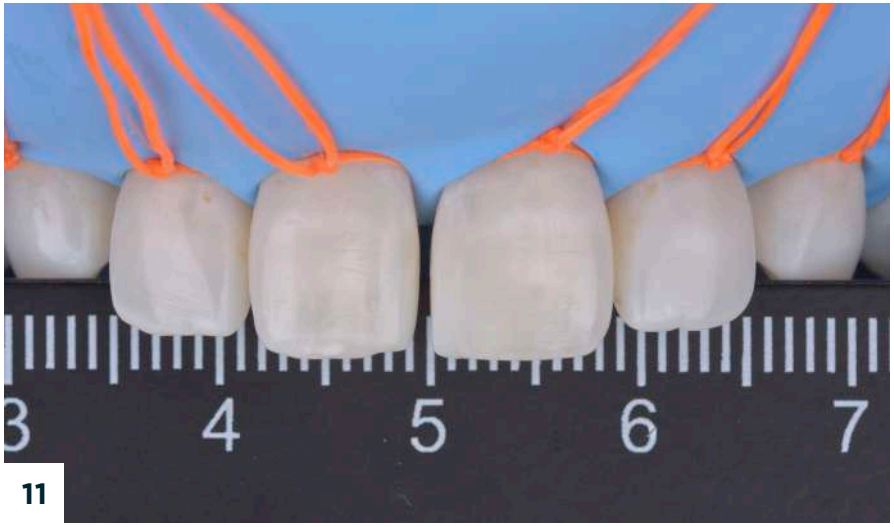


11

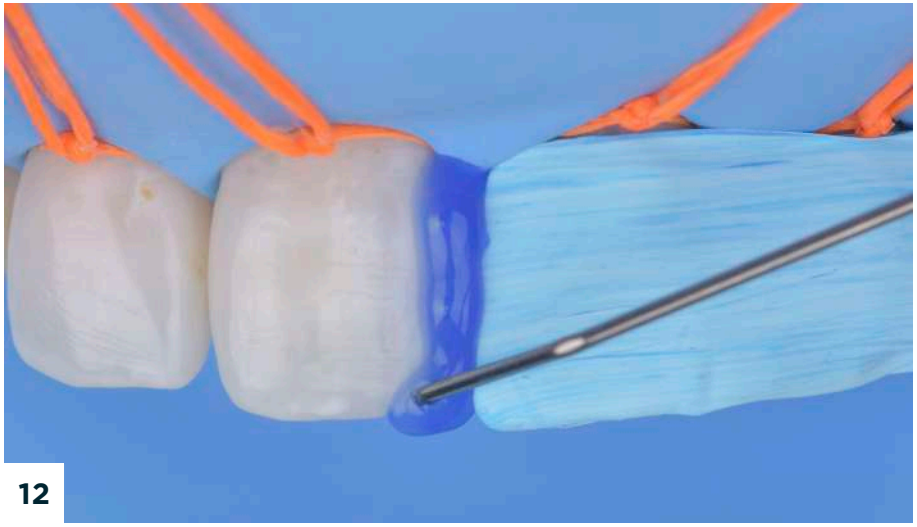

12

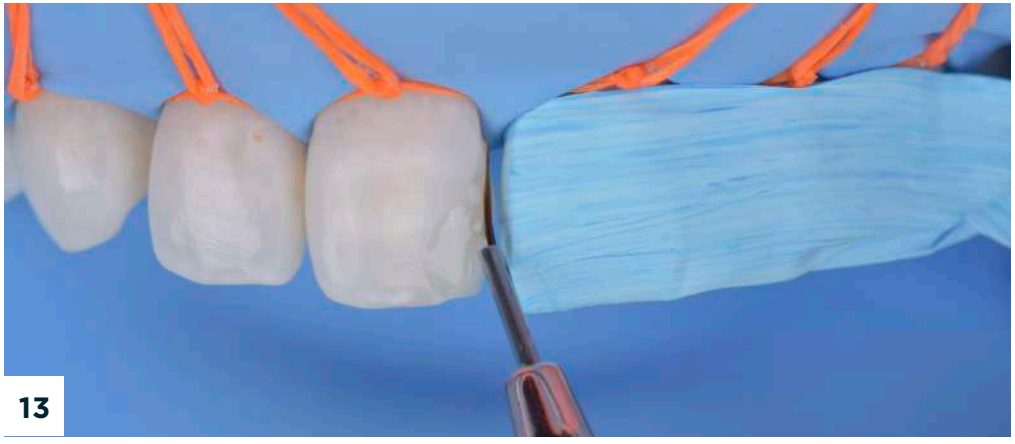

13

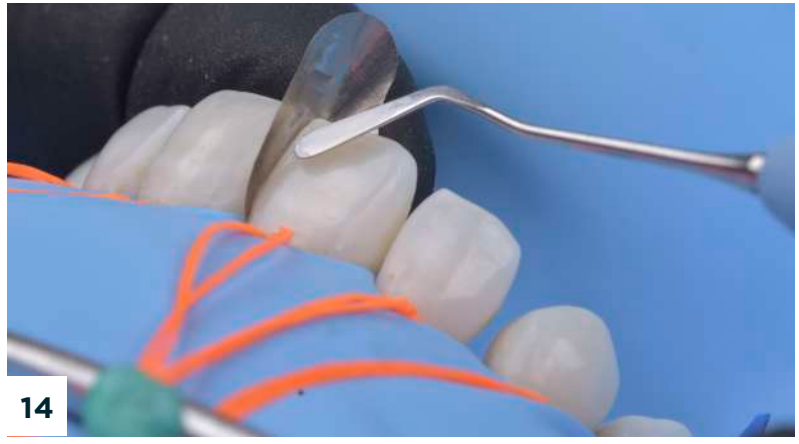


14


11. Verification of size proportion between central incisors.
12. Acid etching on enamel with **Condac 37** on the mesial of tooth 11.
13. Application of **Ambar Universal APS Plus** adhesive.
14. Application of composite increment with **Elóra APS**, in shades B1 on the cervical region and BL on the middle incisal third of tooth 11.




15

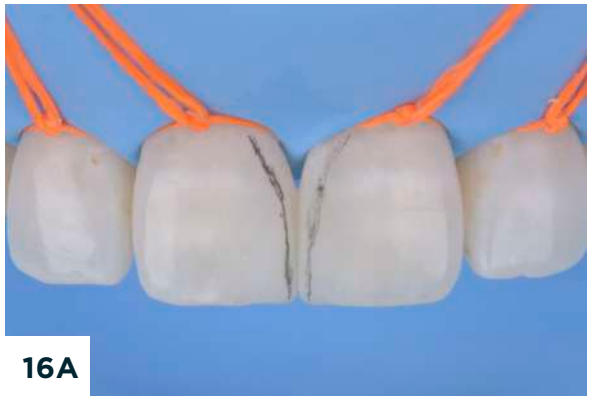

16A

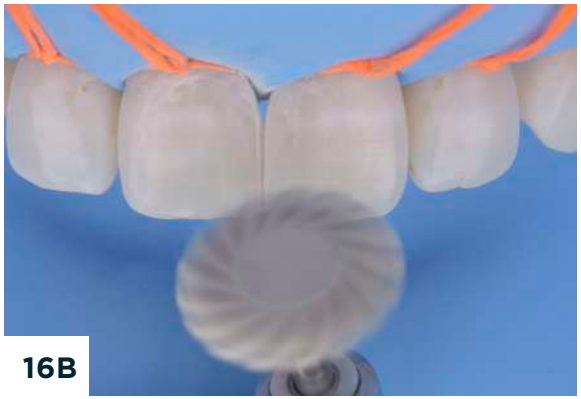


16B

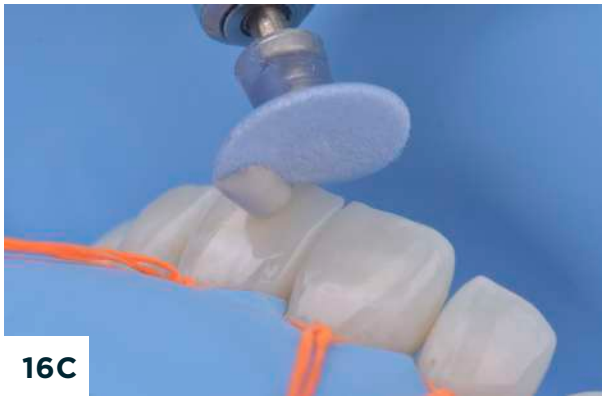


16C


15. Appearence after light-curing.
16A, 16B, and 16C. Finishing and polishing protocol with rubber disc and **Diamond Flex** felt disc, finishing with **Diamond Excel Ultrafine** polishing paste.

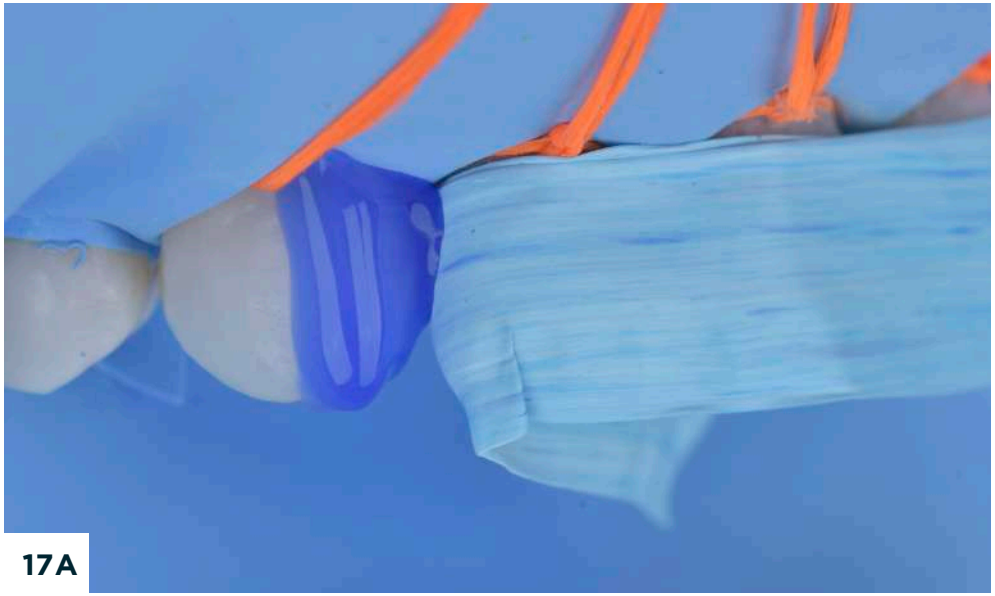

17A

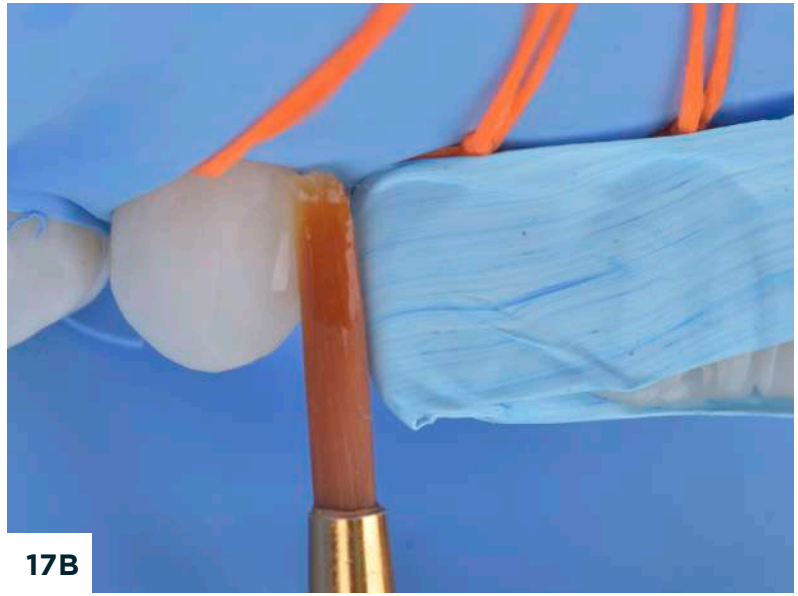

17B

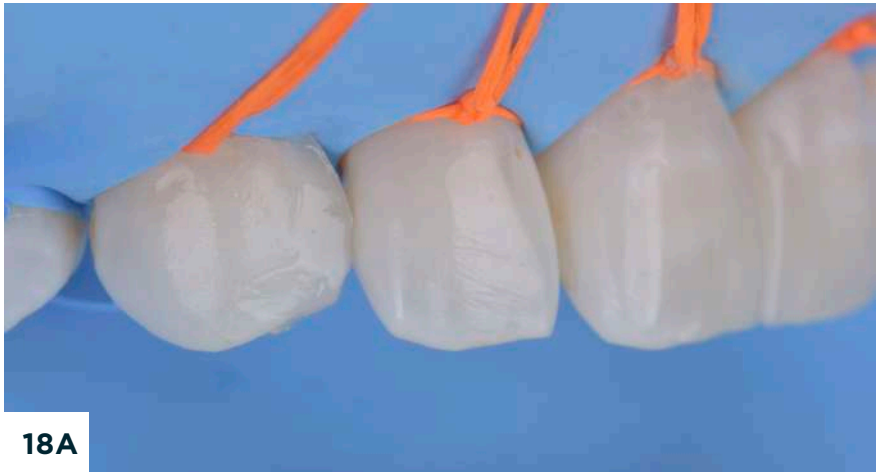


18A



18B

17A and 17B. Acid etching with **Condac 37** on enamel and application of **Ambar Universal APS Plus** adhesive on canines.
18A and 18B. Application of composite increment with **Elóra APS**, in shades B1 and BL, on the mesial of teeth 13 and 23.

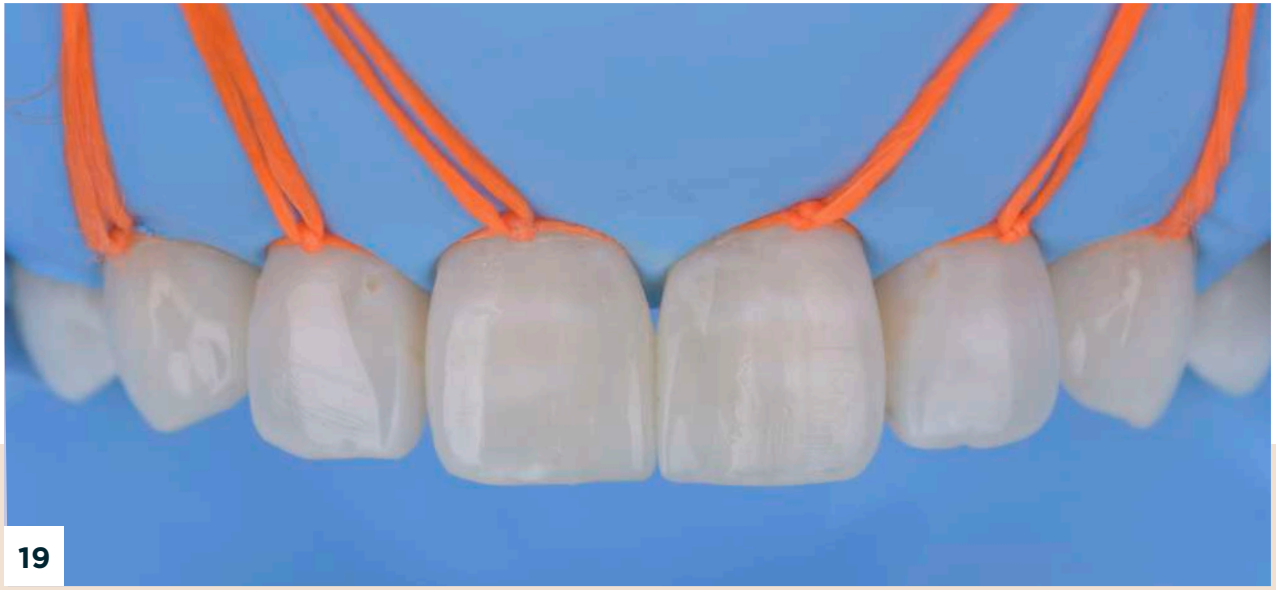

19



20

19. Final appearance after light-curing and finishing and polishing protocol.
20. Final appearance after removal of the absolute isolation.



21

21. Final result of the smile after diastema closure.

Didactic diagram of stratification.



Elóra APS, in shade B1

Elóra APS, in shade BL

A contemporary composite.

One concept, many possibilities..



COMMERCIAL PRESENTATION



Refill:

Syringe with 4 g, in shades:
A1, A2, A3, A3.5, B1, BL, XBL, DB1, DA1, DA3.

Syringe with 2 g in shades:
Trans or WE.

FGM DENTAL GROUP | FGM Elóra by APS

fgmdentalgroup.com

MDSAP MEDICAL DEVICE SINGLE AUDIT PROGRAM

EU-MDR

APPROVED FDA APPROVED

Obelis

BUENAS PRÁCTICAS DE FABRICACIÓN

ISO 13485:2016

1000056629-01